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# Thank you!

Thank you for choosing this BEYOU product. BEYOU products are chosen, tested, and recommended by ELECTRO DEPOT, so you can be sure you are getting a top-quality, easy-to-use product that won't let you down.

We're confident your new device will be a pleasure to use every day!

Welcome to ELECTRO DEPOT.

# Table of Contents

**A**

## **Product overview**

User-friendly Tips  
Features/Specifications  
Performance Index

---

**B**

## **Product usage**

To Install or Replace Battery  
Key Illustration  
LCD Display  
To Set Weight Unit  
Optional Normal Weighing Mode (by passing parameter setting)  
Fat/Hydration/Muscle/Bone/Calorie Measuring Mode  
To recall test record  
Warning Indications  
About Fat  
About Calorie

---

**C**

## **Cleaning and maintenance**

Cleaning and maintenance  
Storage

This product is specially designed to test body fat/hydration/muscle/bone, help keep a close eye on the daily change of your body fat, hydration ,muscle, bone and weight, providing information for your reference. Base on the analysis of your height,weight and age etc., it tells your daily energy requirement, and helps control your weight. Please read this instruction carefully before use.

## Product overview

- Unpack the product. Remove all the labels from the product except the nameplate. Check that it is complete and in a good condition. If the unit is damaged or faulty, do not use it and take it to your retailer or after-sales service centre.
- Keep all packaging out of the reach of children. There is a risk of accidents if children play with packaging materials.

## Intended use

- These scales have been designed to only weigh people with a weight of up to 150 kg. Do not use them for any other purpose!
- Any other use may damage the unit or cause injury.

## Specifications

Model:	981753 - BY-IMPED
Power supply:	2x1.5V AAA batteries
Maximum load:	150kg / 330lb

## User-friendly Tips

### During Measuring

- The device is designed to allow auto step on. Please always ignore the first reading and only start taking effective record from the second time weighing.
- Always remove your shoes and socks , wipe clean your feet before using.
- Step on the weighing platform gently.
- It is advisable to take measurement at the same time of the day.
- Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.

- Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
- The measured data of the following people may have deviation:
  - Children under 10 or adult over 100 (can only use the device under normal weighing mode)
  - Adult over 70
  - Body building athletes or other occupational athletes.
- People with the following conditions should not use the device:
  - Pregnant women
  - Those with symptom of edema
  - Those on dialysis treatment
  - Those who use heart pacemaker or with other implanted medical device.

### Auto-off

If there is no operation in 30 seconds, scale will be auto off.

## Features/Specifications

1. Using high precision strain gauge sensor
2. Multifunction: test body fat%, body hydration%, body muscle%, body bone% and suggest calorie intake requirement.
3. Built-in memory for 12 different users
4. Optional normal weighing mode (bypassing parameter setting): Auto step on weighing and consecutive add-on weighings
5. Auto user recognition function
6. Low battery and overload indication

## Performance Index

Item Index	Weighing Capacity	Fat%	Hydration%	Muscle%	Bone%	Calorie
Index Range		4.0% ~60.0%	27.5% ~66.0%	20.0% ~56.0%	2.0% ~20.0%	
Division	0.1 kg 0.2 lb	0.1%	0.1%	0.1%	0.1%	1 KCAL

## To Install or Replace Battery

The battery type of scale depends on the product structure. Find the battery type of your scale according to the following pictures, then install and replace battery as described.

- 2 x 1.5V (AAA) batteries (included). Remove package foil before use;

Open the battery case cover at the bottom of scale:


- A. Take out used batteries, with the aid of an object with pointed end if needed.
- B. Install new batteries by putting one side of the battery down first and pressing down the other side. (Please observe the polarity)
- C. Reinsert the battery compartment cover until it clicks into place.

## Key Illustration

SET (ON/ SET): Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter overleaf)

- ▲ : Advance the value by one interval each press. Continuous pressing will speed up the advance of value.
- ▼ : Decrease the value by one interval each press. Continuous pressing will speed up the decrease of value.


## LCD Display

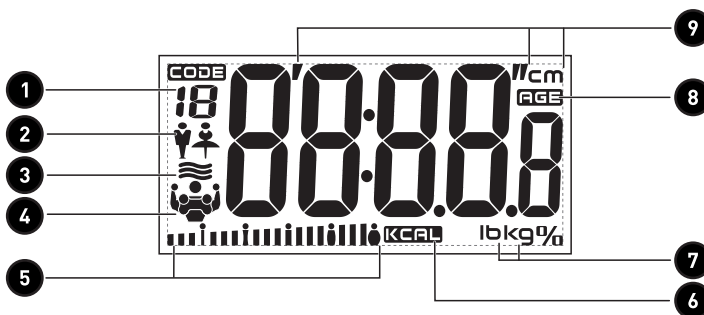1	User memory code	6	Calorie Unit
2	Gender	7	Weight Unit
3	Hydration%	8	Age
4	Fat% / Muscle% / Bone%	9	Height Unit
5	Fat status indication		


## To Set Weight Unit

If your scale is available (or applicable) in kg/lb, you can set weight unit as follow:

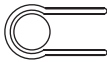At the bottom of scale, you can find a unit conversion button as shown in picture.

Press the unit conversion button when scale is power on, LCD shows current weight unit. Press the button again to convert weight unit.

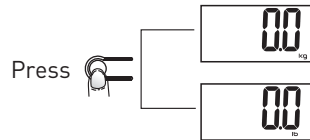
- When setting lb as weight unit, the default height unit is feet/inch;
- When setting kg as weight unit, the default height unit is cm.

**Remark:** The unit system of scale depends on sales destination or the requirement of distributor.

## Optional Normal Weighing Mode (by passing parameter setting)


### Auto step on weighing
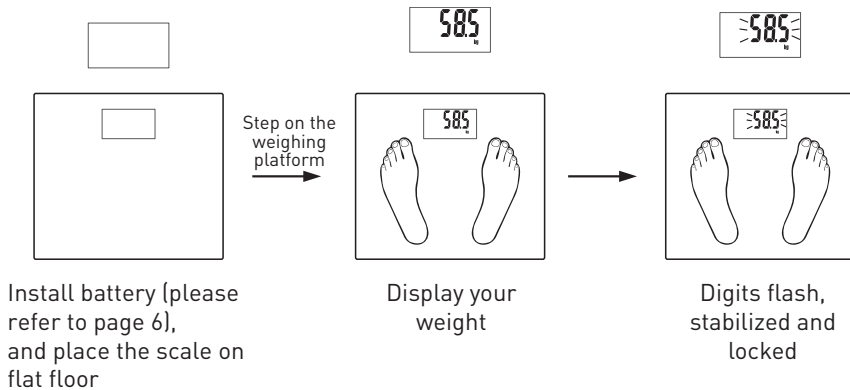
### Consecutive add-on weighing

After digits stabilize, more additional loadings can be conducted. Ensure each add-on loading has to exceed 2kg to get a new reading.

Example: If you want to get a weighing value of your baby, you can weigh yourself first, and then hold your baby together, without stepping off the scale. Then you will get the total weight reading. Subtracting the previous reading of yourself from the new reading, you will get the difference as your baby's weighing value.
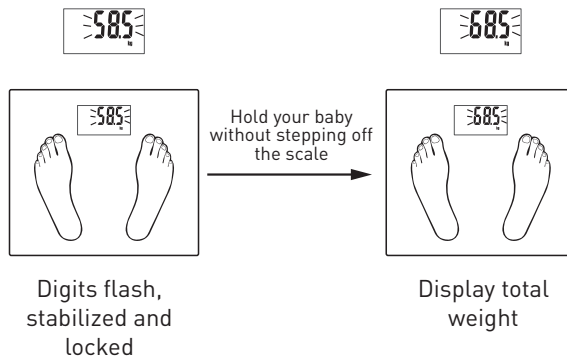
## Fat/Hydration/Muscle/Bone/Calorie Measuring Mode

### Setting personal parameters

- At the first time when you install battery and switch on, all default parameter are set as follow:


Memory	Gender	Height	Age
1	Male	165cm / 5'05.0"	25

- Parameters can be set in the following range:

Memory	Gender	Height	Age
1~12	Male / Female	100 ~ 250cm 3'03.5" ~ 8'02.5"	10 ~ 100

- Parameter setting

Example: the 3rd user, female, 168cm, 28 years old, can set her parameters as follow:

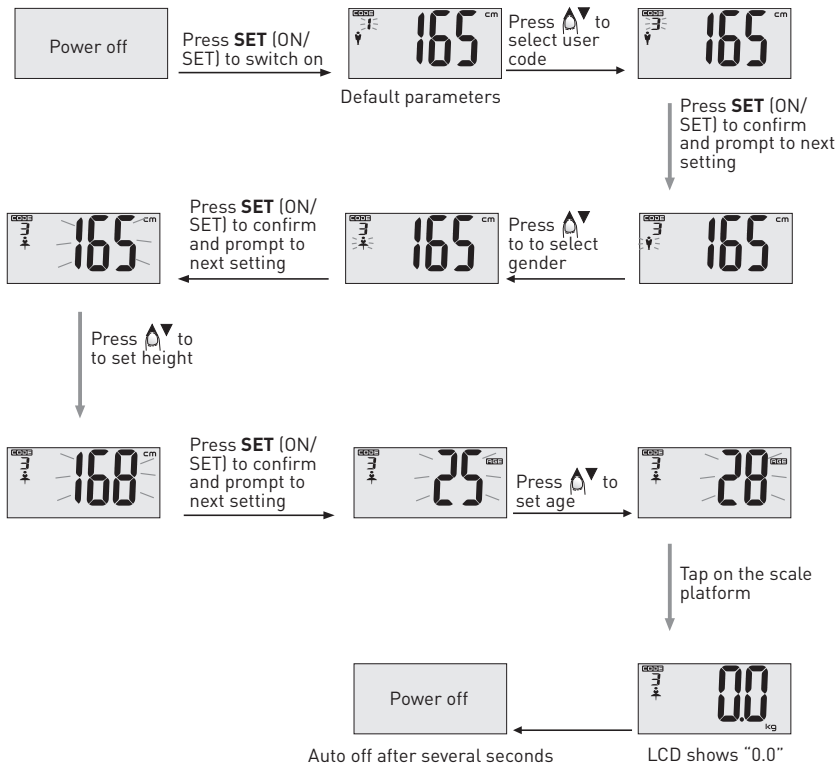
**NOTE:**

- During parameter setting, if idle for 6 seconds, the parameters already set will be memorized automatically and LCD shows "0", meaning body fat measuring will be started.
- If LCD prompt to show "0" during parameter setting, press SET (ON/SET) for 3 seconds to return setting mode, then you can set user code, gender, height and age.


**To Start Measuring** (Make sure your feet are touching with metal electrodes on the scale platform, or body fat can not be measured.)

To ensure accuracy, please remove your shoes & socks before stepping on. Example: the 3rd user, female, 168cm, 28 years old can measure as follow:
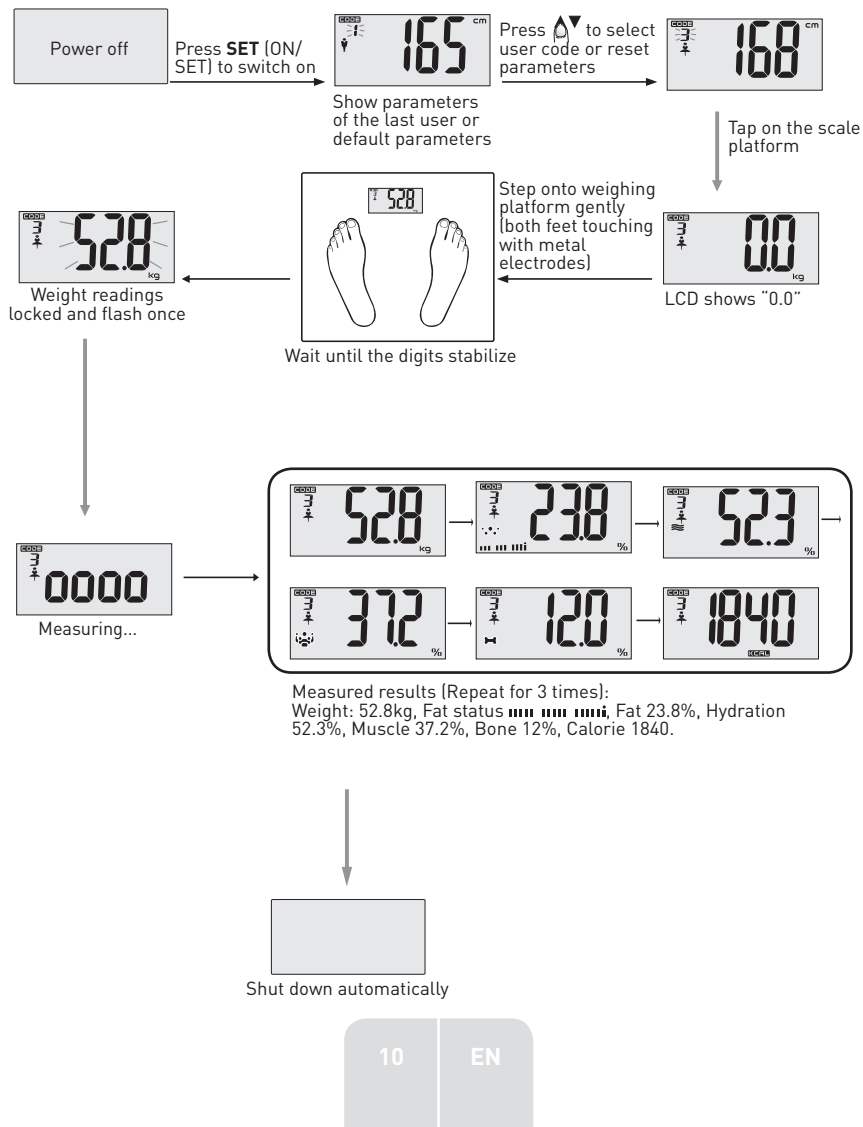
### Auto User Recognition


When scale is power off or standby, please step on the scale to get a weighing value. LCD will lock the reading in one second. If the scale can search similar weighing value (value difference below  $\pm 2\text{kg}$ ), LCD will flash weighing value once and show user code. The scale will automatically recognize current measuring user after flashing user code once. Then scale will enter into measuring other data, like body fat, hydration, etc.
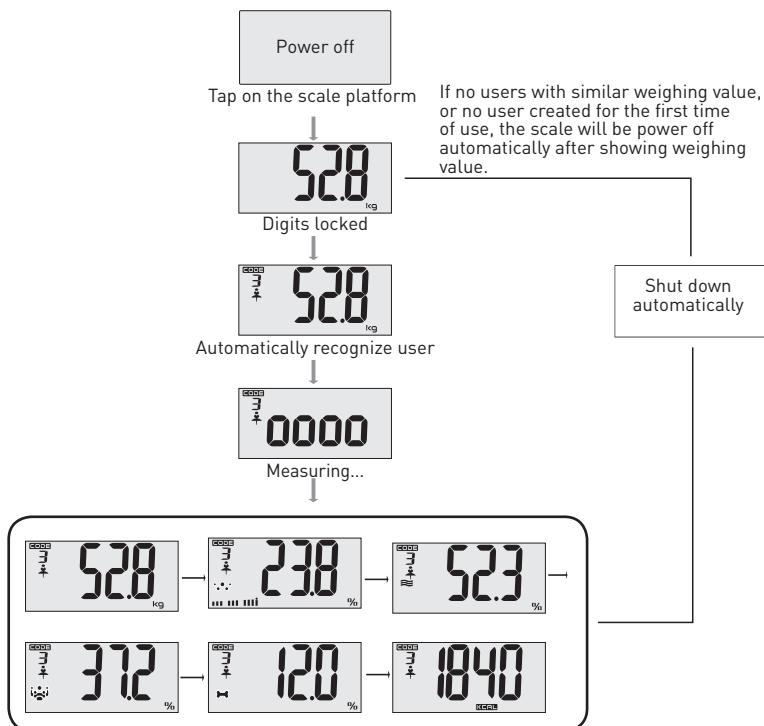
#### Remark:

- After searching users with similar weighing value and flashing the weighing value once, the scale can not enter into body fat measuring if users do not touch metal electrodes effectively. The scale will only show users' gender and current weighing value, till automatically shut down.
- If the scale find two or more similar weighing value, it will choose the user with closest weighing value. If the scale find two same weighing values, it will choose the user code in front.

## To recall test record
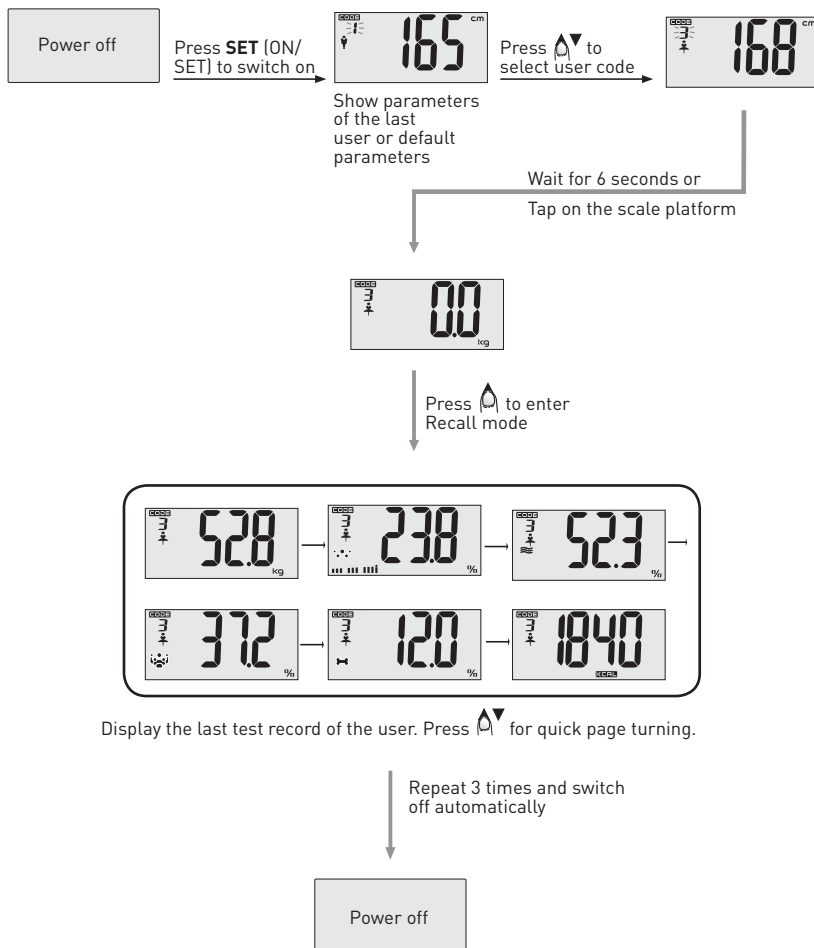
## Warning Indications

### Low Battery Indication

The battery power is running low, please replace with a new battery.

### Over-load Indication

The weighing subject on the platform exceeds the maximum capacity of scale. Please step off to avoid damage.

### Low Fat% Indication

The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

### High Fat% Indication

The Fat% is too high. Please watch your diet and do more exercise.

### Re-test

Error exists, please re-test to get the correct result.

## About Fat

### 1. About Fat%

- Fat% is an index % of fat content in human body.
- The fat content in human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may invade. Here goes an old saying, "Obesity is hotbed of illness".

## 2. Working Principles

As electrical signal of certain frequency transmits through human body, the impedance of fat is stronger than that of muscle or other human tissues.

Through the use of biological resistance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signal transmitting through human body.

## 3. Fat/Hydration Content Fitness Assessment Chart

After your weighing, you can use the following Assessment Chart for reference:

Age	Fat (F)	Hydration (F)	Fat (M)	Hydration (M)	Status	Key Line Indication
≤30	4.0-16.0%	66.0-57.8%	4.0-11.0%	66.0-61.2%	Too lean	
	16.1-20.5%	57.7-54.7%	11.1-15.5%	61.1-58.1%	Lean	
	20.6-25.0%	54.6-51.6%	15.6-20.0%	58.0-55.0%	Normal	
	25.1-30.5%	51.5-47.8%	20.1-24.5%	54.9-51.9%	Fat	
	30.6-60.0%	47.7-27.5%	24.6-60.0%	51.8-27.5%	Too fat	
>30	4.0-20.0%	66.0-55.0%	4.0-15.0%	66.0-58.4%	Too lean	
	20.1-25.0%	54.9-51.6%	15.1-19.5%	58.3-55.3%	Lean	
	25.1-30.0%	51.5-48.1%	19.6-24.0%	55.2-52.3%	Normal	
	30.1-35.0%	48.0-44.7%	24.1-28.5%	52.2-49.2%	Fat	
	35.1-60.0%	44.6-27.5%	28.6-60.0%	49.1-27.5%	Too fat	

The above information is only for reference.

## About Calorie

### 1. Calorie & Energy Consumption

Calorie: Energy Unit

The data reading indicates the energy you need to sustain the daily metabolism and ordinary activities required for your weight, height, age and sex.

### 2. Calorie & Losing Weight

Control the calorie intake combined with proper exercise is an effective way to lose weight. This is a very simple theory. As the daily energy intake can not support the energy consumption of the human body, the saccharide and fat stored inside will decompose and provide energy to human body.

That leads to weight loss.

## Cleaning and maintenance

- The platform will be slippery when it is wet. Keep it dry!
- Keep standing still during the course of weighing.
- DO NOT strike, shock or drop the scale.
- Treat the scale with care as it is a precise instrument.
- Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/abrasive cleaning agents.
- Keep the scale in cool place with dry air.
- Always keep your scale in horizontal position.
- If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
- If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault. If you can not resolve the problem, please contact your dealer for consultation.

## Storage

- Store the unit in a fresh, dry and dust-free place, out of the reach of children and pets.
- Please retain the packaging and use it to store your unit when it is not in use for an extended period of time.

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